Horizon™

HORIZON 627 LUMBAR, 631 LSO LoPro, 637 LSO

Aspen
MEDICAL PRODUCTS

ASPEN MEDICAL PRODUCTS

6481 OAK CANYON, IRVINE, CA 92618-5202, USA
800-295-2776 · 949-681-0200 · aspenmp.com

ENGLISH

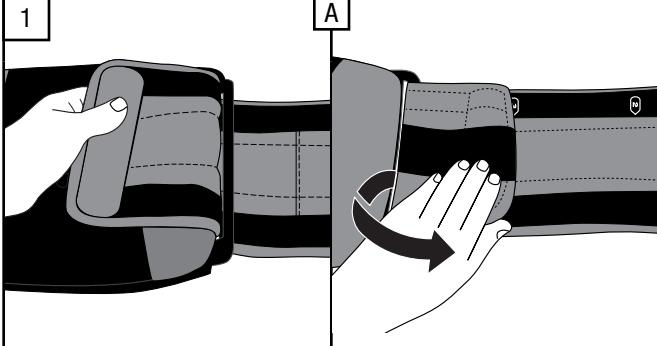
ESPAÑOL

DEUTSCH

ITALIANO

FRANÇAIS

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
| SIZING |
|---------------------------------------|
| CIRCUMFERENCE |
| U.S. MEN'S PANT SIZE RANGES |
| U.S. WOMEN'S WOMEN'S PANT SIZE RANGES |

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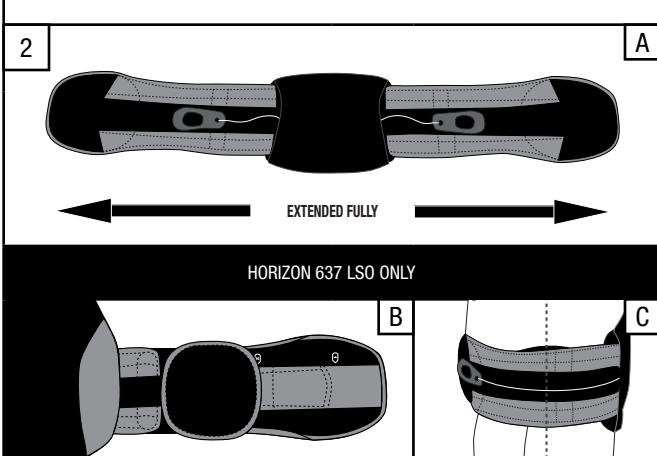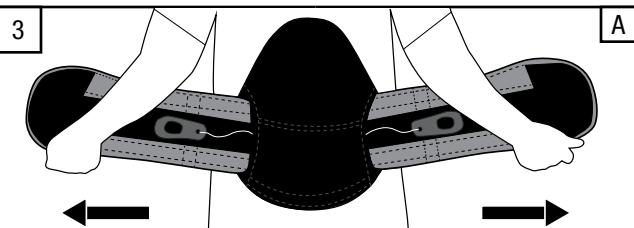
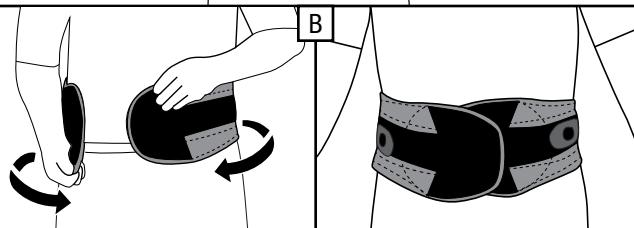
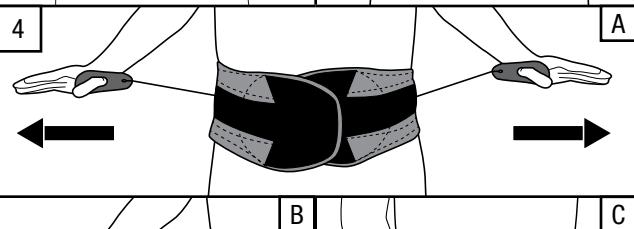
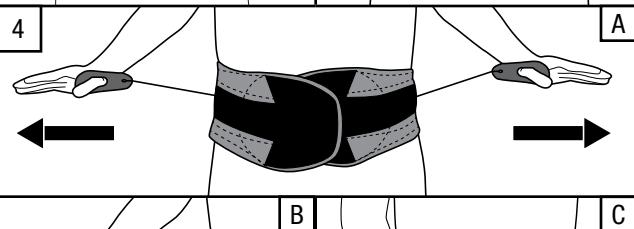
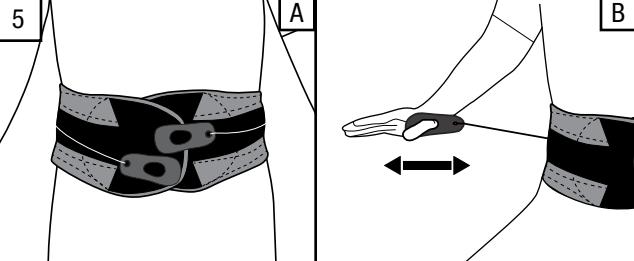
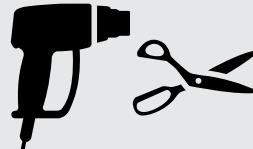
THIS PRODUCT MAY BE SUBSTANTIALLY MODIFIED BY AN INDIVIDUAL WITH EXPERTISE

For more information, please visit
www.aspenmp.com/education/product-modifications

LIMITED WARRANTY

Aspen Medical Products, Irvine, CA 92618, warrants to the user who originally purchases this product that it is free from defects in material and workmanship. The sole obligation of Aspen Medical Products in the event of breach of warranty shall be to repair or replace the defective product or part(s).

Aspen Medical Products shall have no obligation under this limited warranty in the event:

- (a) The product was not purchased from Aspen Medical Products or through its authorized channels of distribution;
- (b) The product is altered;
- (c) Any parts not supplied by Aspen Medical Products are inserted into the product; or
- (d) The product is not used in accordance with the Aspen Medical Products Instructions for Use.

THE FOREGOING IS THE SOLE AND EXCLUSIVE REMEDY FOR FAILURE IN SERVICE OF, OR DEFECT IN, THE PRODUCT. ASPEN MEDICAL PRODUCTS SHALL NOT BE LIABLE UNDER THIS OR ANY IMPLIED WARRANTY FOR ANY DIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES. THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING THE WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE OR USE, AND ALL OBLIGATIONS OR LIABILITIES ON THE PART OF ASPEN MEDICAL PRODUCTS FOR DAMAGES ARISING OUT OF OR IN CONNECTION WITH THE USE OF THE PRODUCTS, WHICH ARE HEREBY DISCLAIMED AND EXCLUDED BY ASPEN MEDICAL PRODUCTS.

READ INSTRUCTIONS BEFORE USE. Proper training in the use of this device should take place before it is applied. These directions are guidelines only and are not offered as medical recommendations.

THIS PRODUCT IS INTENDED FOR APPLICATION BY LICENSED HEALTH CARE PRACTITIONERS AS DIRECTED BY A PHYSICIAN OR OTHER QUALIFIED MEDICAL AUTHORITY. PERMANENT INJURY MAY RESULT FROM USE OTHER THAN AS DIRECTED BY A PHYSICIAN OR OTHER MEDICAL AUTHORITY. ASPEN MEDICAL PRODUCTS SPECIFICALLY DISCLAIMS ANY OBLIGATION OR LIABILITY FOR INJURY OR DAMAGE BY REASON OF ANY APPLICATION OF THIS PRODUCT. PLEASE SEE THE LIMITED WARRANTY FOR FURTHER INFORMATION IN THIS REGARD.

This warranty gives you specific legal rights and you may have other rights which vary from state to state.
 Warning: To ensure proper use, please review all material.
 A listing of authorized distributors in your area is available upon request.